

Personal -86°C ULT Freezer



DW-86L100J

Typical Installation and Application

Suitable for clinical, medical, scientific research, quarantine and other departments to store items under low temperature conditions. Applicable for universities, hospitals, disease prevention and control centres, blood stations, scientific research institutes, electronics and chemical enterprise laboratories and biomedical engineering research institutes.

For storage of biological products and biological samples such as red and white blood cells, viruses, bone and bacteria. Also used for electronic devices and other materials used for cryogenic tests.



Energy Efficient, Safe and Reliable

High efficiency HC refrigeration system, optimised for energy efficiency delivering a power consumption figure of just 5.5 kW·h/24 h



Personal ULT Storage

810mm cabinet height makes it easy to place on or under bench, saving storage space. Stackable design.



Ergonomic Design

Ergonomic handle design ensures easy one-hand door opening.



Low Noise

Optimized noise reduction cabinet and system design, emits sound level of only 46.8 dB.

Ihr Fachhändler:
www.axon-labortechnik.de

Axon
Labortechnik

Axon Labortechnik GmbH
Im Braumenstück 16 - D-67659 Kaiserslautern
Verkaufsbüro Tel. 06223/95 41 190 Fax: 06223/95 40 714
info@axon-labortechnik.de



Haier Biomedical
International



Haier Biomedical
International



@haiermediclnt



Haier Biomedical
International




Haier Biomedical
International

Product Features



Haier Biomedical -86°C Personal ULT Freezer Energy-efficient HC refrigeration design. Secure and reliable. Intelligent Control.

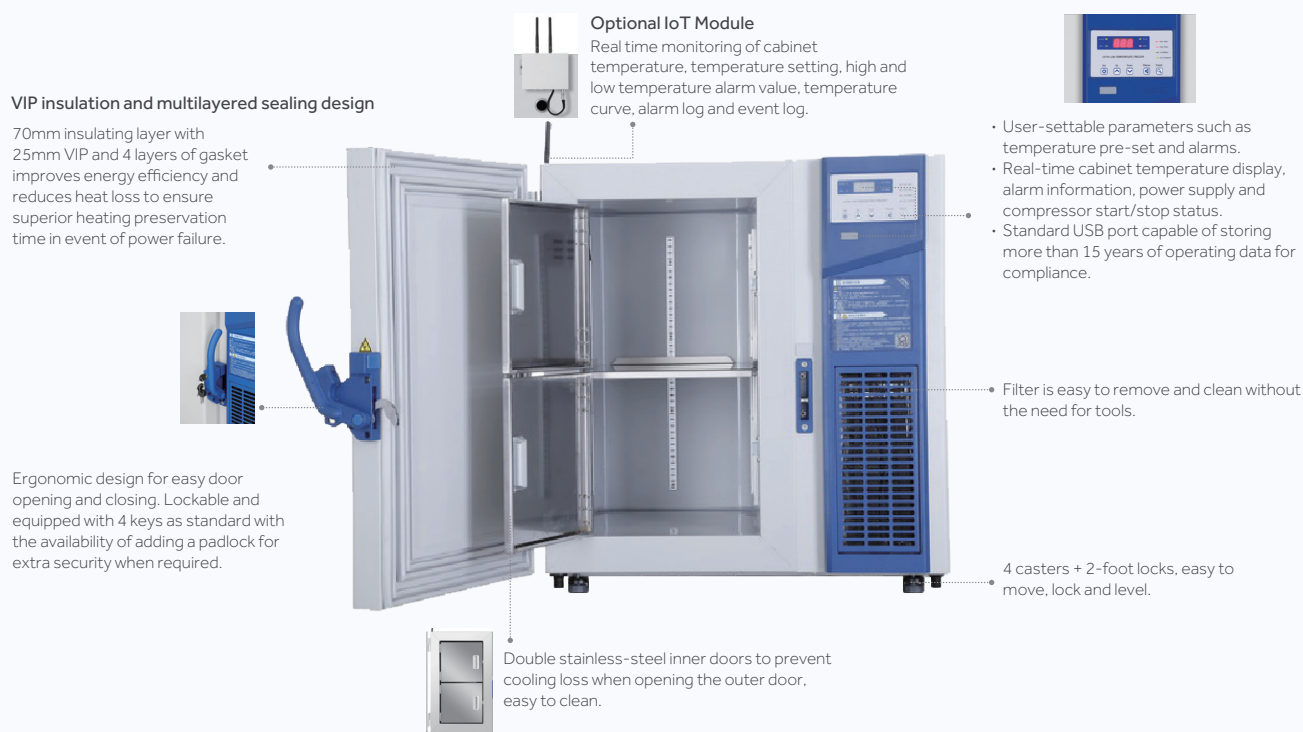

Microprocessor control system

- Microprocessor control with LED temperature display. The adjustable precision is 1 °C. And the temperature in the unit is adjustable from -40 °C to -86 °C.
- Cabinet temperature/voltage/ambient temperature checking are available.
- Multiple fault alarms: High/low temperature, high/low voltage, high ambient, sensor error, door ajar, power failure, low battery, remote alarm, and hot condenser alarm.
- Two types of alarms: Audible buzzer alarm and flashing light alarm.
- Battery backup alarm function operates continuously for >24hrs in the event of a power outage.
- Standard configuration: RS485 port and USB interface.
- Standard 5V power supply available for test equipment.
- Optional IoT module.



Superior thermal insulation performance

70mm super thick insulation layer design, aviation vacuum insulation material VIP, thickness of 25 mm or more, 4 layers of silicone seal, superior thermal insulation and energy saving effect.



Porthole

Portholes as standard, allows for independent testing of cabinet temperature.



Security lock

Standard door lock and padlock to ensure sample security and prevent unauthorised access.



USB data storage




Capable of storing more than 15 years of data.



Typical Performance Characteristics in 25 °C Ambient: DW-86L100J



— Ambient Temp.
— Cabinet Temp.

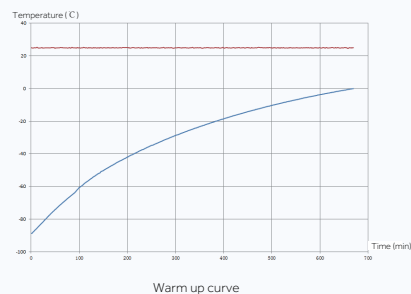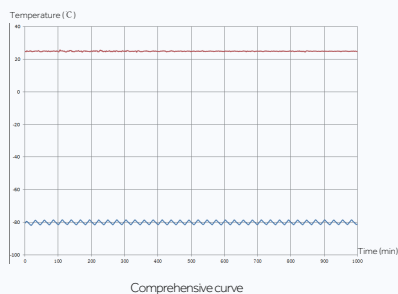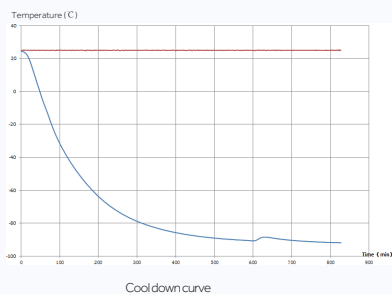

**Cooling Time(ambient temperature
reducing to -80 °C)**

310min



Temperature Uniformity

±4°C

**Time for inside temperature to rise back
from -80 °C to -50 °C naturally at ambient 25 °C**

120min

Product Dimension Drawing

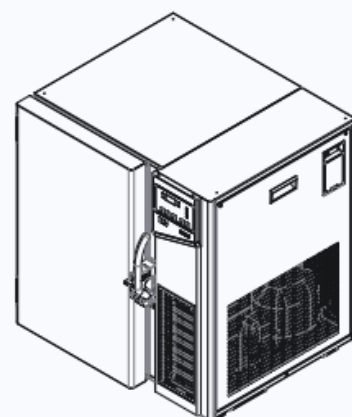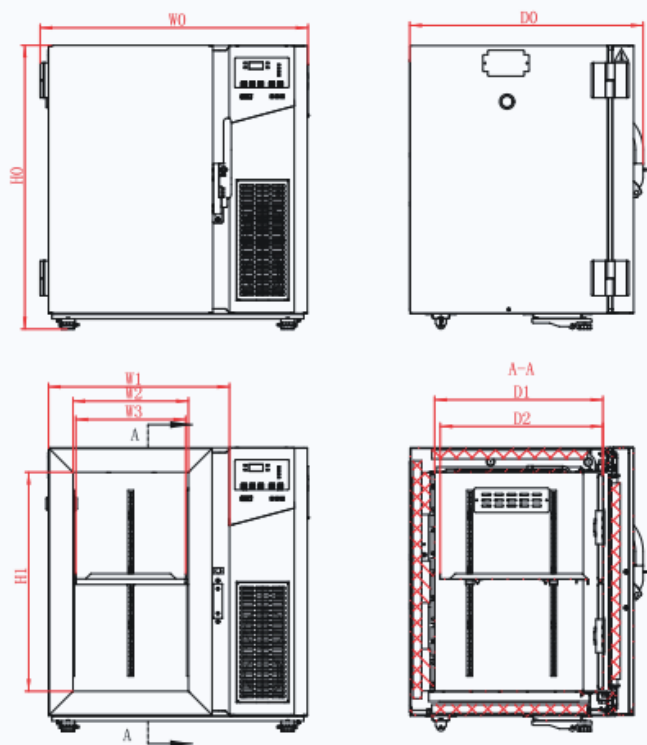

KEY

W0=TOTAL WIDTH
W1=FREEZER WIDTH
W2=INNER CABINET WIDTH
W3=INNER CABINET SHELF WIDTH
D0=TOTAL DEPTH
D1=INNER CABINET DEPTH
D2=INNER CABINET SHELF DEPTH
H0=TOTAL HEIGHT
H1=INNER CABINET HEIGHT

CODE MODEL	W0	W1	W2	W3	D0	D1	D2	H0	H1
DW-86L100J	770mm	520mm	330mm	314mm	660mm	481mm	465mm	810mm	630mm



Model		DW-86L100J	
Technical Data	Cabinet Type	Upright	
	Climate Class	N	
	Cooling Type	Direct cooling	
	Defrost Mode	Manual	
	Refrigerant	HC	
	Noise (dB (A))	46.8	
Performance	Cooling Performance (°C)	-86	
	Temp Range (°C)	-40~-86	
Control	Controller	Microprocessor	
	Display	LED	
Electrical Data	Power Supply (V/Hz)	220~240/50	120/60
	Power (W)	680	680
	Electrical Current (A)	3	6.5
	Capacity (L/Cu.Ft)	100/3.5	
Construction	Net/Gross Weight (approx)	(kg)	108/132
		(lbs)	238/291
	Interior Dimension (W*D*H)	(mm)	330*481*630
		(in)	13*19*25
	Exterior Dimension (W*D*H)	(mm)	770*660*810
		(in)	30*26*32
	Packing Dimension (W*D*H)	(mm)	830*710*970
		(in)	32*28*38.5
Loading Quantities	Container Load (20'/40'/40'H)	44/88/88	
Alarms	Remote Alarm	Y	
	High/Low Temperature	Y	
	Hot Condenser	Y	
	Power Failure	Y	
	High/Low Voltage	Y	
	Sensor Error	Y	
	Low Battery	Y	
	High Ambient Temperature	Y	
	Door Ajar	Y	
Accessories	Caster	Y	
	Foot	Y	
	Porthole	Y/1	
	Shelves/ Inner Doors	1/2	
	USB Interface	Y	
	5V Power Supply Port	Y	
	RS485 Port	Y	
Certification	CE/UL	Y/Y	

Product appearance and specifications are subject to change without notice