

Slide Spinner

The *Slide Spinner* is the perfect centrifuge for drying microarray slides before scanning.


Drying of microarray slides is critical for producing quality imaging and screening results. The *Slide Spinner* dries 2 slides in as little as 10 seconds. No set up is necessary - simply place the slides into the cassettes and close the lid of the centrifuge to start. The unit quickly accelerates to 4,800 rpm, the optimal speed for drying the slides without splashing. Upon opening the lid, the rotor comes to a stop.

All liquid from the slides is captured in the cassettes. After centrifugation, the cassettes are easily removed and the liquid poured off. The cassettes can then be rinsed and placed back in the centrifuge.



CE

- Quickly dries microarray slides
- Capacity for two slides
- Removable cassettes catch all liquid
- No splashing
- Starts and stops in seconds

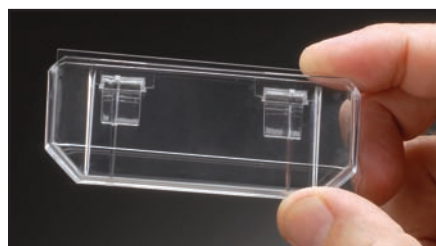

SPECIFICATIONS

Speed	4,800 rpm
Maximum capacity	2 standard slides
Dimensions (W x D x H)	5.9 x 5.9 x 4.6 in/15 x 15 x 11.7 cm
Weight	1 lb/0.45 kg
Electrical	230V~, 50 Hz or 120V~, 60 Hz

ORDERING INFORMATION

Catalog No.	Description
C1303-T-230V*	Slide Spinner with two slide drying cassettes, 230V
C1303-SC	Extra cassettes for Slide Spinner centrifuge, pack of 2

*When ordering 230V units, please specify cord type by adding -EU or -UK to the end of the catalog number. 120V units are available. When ordering, please specify by removing -230V from catalog number.

Ihr Fachhändler:
www.axon-labortechnik.de

Axon
Labortechnik

Labnet